# **INFINITY2-3**




### 3.3 Megapixel CCD USB 2.0 Camera

High Resolution CCD Color
Microscopy Camera for Qualitative Image Archiving

### **Outline**

INFINITY2-3 digital CCD color camera offers excellent sensitivity and high resolution. This is an ideal camera for brightfield, darkfield, and DIC imaging applications where higher resolution archiving and a wide dynamic range are routinely required. With 2080x1536 resolution and on-board processing, the INFINITY2-3 delivers outstanding image quality for a wide variety of scientific applications.

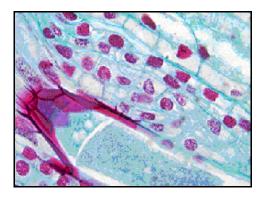

Live video preview provides for real-time focus while auto exposure and auto white balance efficiently capture your optimal image. An intuitive user application provides camera controls while full integration to popular third-party imaging applications is available through our TWAIN drivers.

The USB 2.0 digital interface ensures a simple plug-and-play installation – and one standard cable minimizes camera clutter. No framegrabber required.

Every INFINITY camera includes INFINITY ANALYZE software for advanced camera control, image processing, measuring and annotation.

The INFINITY2-3 camera is supported by an experienced team of technical support and imaging experts. We understand your imaging needs and are here to help you get the most out of your camera.

### **Performance Features**

- ☐ The high-speed USB 2.0 interface eliminates a framegrabber and facilitates ease of installation on both laptop and desktop computers
- ☐ The low noise characteristic of the INFINITY2-3, 3.3 megapixel CCD image sensor results in crisp color quality for the most demanding brightfield and darkfield microscopy applications including clinical pathology and cytology, histology, life science and geology
- □ Full color sub-windowing allows for rapid focus and scanning of samples: 5 fps at full 2080X1536 resolution
- ☐ Select 8 and 12-bit pixel data modes
- ☐ The RGB data captured through each pixel contains 36-bits of color image information resulting in 4096 intensity values
- Camera control through an intuitive user TWAIN interface results in rapid image capture archiving and documentation for high throughput applications, demanding research environments and teaching facilities.
- The INFINITY2-3 has a compact design equipped with a C-mount, facilitating installation on all microscope configurations including upright, inverted and stereo
- □ INFINITY 2-3 cameras are software compatible with Windows XP, Windows Vista, Windows 7, Mac OSX, 32 and 64-bit

## **Specifications**

### **INFINITY2-3**

#### Camera Sensor

Image Sensor 1/1.8" Sony ICX262 interlaced scan 3.3 megapixel

color interline CCD sensor

Effective Pixels 2080 x 1536, 3.45 µm square pixels

Frame Rate 5 fps at full resolution, higher fps with binning and

ROI

Digital Output 8 and 12-bit
Read Noise 12 e- rms
Readout Frequency 22.5 MHZ

Dimensions (W x H x D) 2.25 x 3.85 x 1.56 inch (enclosed)



Mass ~300g

Power Requirement USB bus power, or external 5 V DC – 500 mA

Power Consumption ~2.5 W

Operating Temperature 0 to +50 ℃

Operating Humidity 5 to 95 %, Non-condensing

Integration Time 1/1000 to 16 sec.

Shutter Global

ROI User selectable through mouse click and drag

Auto Exposure

Automatic / Manual

White Balance

Gain

Automatic / Manual

Programmable / 1to 10x

Interface Connector Standard USB 2.0 high-speed interface

Dimensions (L x W x H) 3.85 x 2.00 x 2.75 inch

Lens Mount C-mount lens adapter

Binning Options 2x2, 3x3, 4x4

#### Recommended PC Specs:


- Pentium 4, 1.3 GHz or higher
- 512 MB RAM
- 60 MB hard drive free space or more
- USB 2.0 Port
- Windows 2000 or XP

#### Minimum PC Specs:


- 600 MHz Processor
- 256 MB of SDRAM
- 20 MB hard drive free space
- USB 2.0 Port
- Windows 98 or ME


#### Product Includes:

- INFINITY2-3 digital camera for USB 2.0
- CD-ROM with INFINITY user application software
- TWAIN driver
- Documentation
- USB 2.0 cable

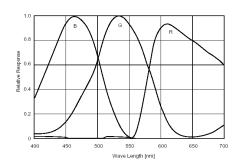

Color Response Curve

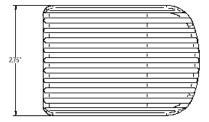

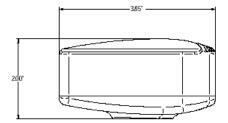

#### Ordering Information

INFINITY 2-3C - Color Camera

Full customization available to meet your exact needs!